


ShortCut™
Split with Ease
Replace with Confidence

ShortCut™ Specifications

ShortCut™ is a groundbreaking solution that mechanically splits pre-existing valve leaflets prior to TAVR, aiming to reduce the risk of coronary obstruction and enable a safer valve-in-valve procedure

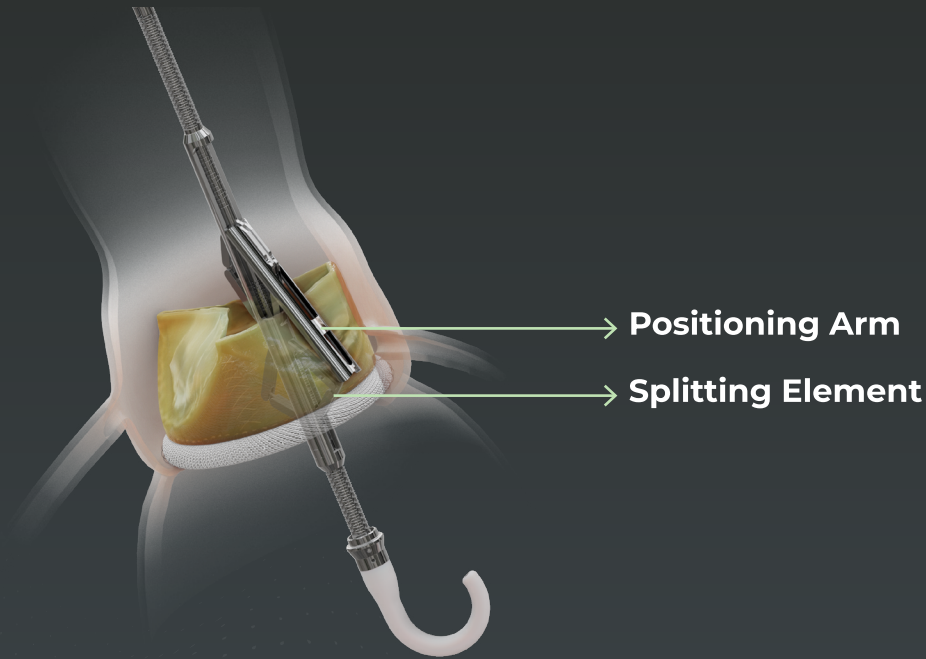


Transfemoral 16 Fr delivery system

Working length: 1150 +/- 15 mm



Compatible with 0.035" guidewire

Ordering Information

Barcode	Product Name	Product Code	Product Description
 (01)7290020101242	ShortCut™ Device	Z-A-20.003-US	ShortCut™ Transcatheter Leaflet Splitting Device

The ShortCut is designed for single use only. Do not re-use or re-sterilize this product. Refer to Instructions for Use for a complete list of warnings, precautions, and contraindications.

ShortCut™ is a medical device cleared by the U.S. Food and Drug Administration for splitting bioprosthetic valve leaflets for patients at risk for coronary obstruction.

PRS-0002677

Pi-Cardia, Inc
1 Broadway St
Cambridge, MA 02142
T (888) 705-5366
E customerservice@pi-cardia.net
www.pi-cardia.net

